

Welch Allyn Stethoscope Specifications

PRODUCT DESCRIPTION	CARDIOLOGY GRADE		PROFESSIONAL GRADE			DISPOSABLE
	HARVEY DLX (28")	HARVEY ELITE DOUBLE HEAD (28") (22")	PROFESSIONAL ADULT	PROFESSIONAL PEDIATRIC	LIGHTWEIGHT	UNISCOPE
PRODUCT NUMBERS	5079-321 Triple, Black 5079-322 Triple, Burgundy 5079-323 Triple, Navy 5079-324 Triple, Green 5079-325 Double, Black 5079-326 Double, Burgundy 5079-327 Double, Navy 5079-328 Double, Green	5079-122 Elite, Black (22") 5079-125 Elite, Black 5079-270 Elite, Burgundy 5079-271 Elite, Navy 5079-284 Elite, Green 5079-125P Elite, Ped Blk	5079-135 Black 5079-137 Navy 5079-289 Blue 5079-139 Burgundy 5079-285 Forest	5079-145 Black 5079-147 Navy 5079-149 Burgundy 5079-287 Forest	5079-73 Velvet Black 5079-74 Poppy Red 5079-75 Dawn Gray 5079-76 Misty Blue	17461 Adult, Yellow 17462P Pediatric, Blue
PRODUCT TYPE	Double/Triple Head	Double Head	Double Head	Double Head	Double Head	Single Head
HEAD DIAMETER INCHES (mm)	1.9 (48.3)	1.9 (48.3)	1.75 (45)	1.2 (30)	1.9 (48.3)	1.9 (49)
OVERALL LENGTH INCHES (mm)	28 (711)	28 (711)	28 (711)	28 (711)	32 (812.5)	29.5 (750) 27 (686)
TUBING LENGTH INCHES (mm)	17 (431.8)	17 (431.8)	21.5 (546.1)	21.5 (546.1)	24 (609.6)	19 (483)
LATEX FREE (Y/N)	Yes	Yes	Yes	Yes	Yes	Yes
CHESTPIECE MATERIAL	Brass	Stainless Steel	Stainless Steel	Stainless Steel	Chrome Plated Zinc	ABS
DIAPHRAGM MATERIAL	Polycarbonate	Polycarbonate	Clear Polycarbonate	Iten Resiten G-10/FR4T Glass Epoxy Sheet	PVC	Polypropylene
DIAPHRAGM STYLE	Flat or Corrugated	Flat or Corrugated	Flat	Flat	Flat	Flat
Adjustable Binaural (Y/N)	Yes	Yes	Yes	Yes	Yes	No
EARTIPS (SOFT/HARD)	Soft/Hard	Soft/Hard	Soft/Hard	Soft/Hard	Hard	Hard
TYPE OF TUBING (SINGLE/DUAL LUMEN)	Dual Lumen	Dual Lumen	Single Lumen	Single Lumen	Single Lumen	Single Lumen
ANTI-CHILL RING (Y/N)	Yes (Brass)	Yes	Yes	Yes	Yes	No
SHIPS ASSEMBLED (Y/N)	Yes	Yes	Yes	Yes	No	Yes
PACKAGING	Printed Box	Printed Box	Printed Box	Printed Box	Plain White Box 1-piece	Plastic Bag
DESCRIPTION	Designed to meet the demanding needs of cardiologists, this stethoscope has fully adjustable binaurals and triple leaf spring for comfort. The head is made of trumpet brass because of its ability to transmit sounds accurately. The dual-lumen tubing is latex-free.		Meets all diagnostic needs, from detecting Korotkoff sounds of blood pressure to analyzing subtle heart murmurs. Stainless-steel chestpiece intensifies body sounds transmitted through tissues. Comfort sealing eartips depress and conform to the shape of ear canals for soft yet airtight fit. Latex-free.		Quality design for blood pressure measurement and general auscultation. Lightweight and flexible single-channel tubing provides comfortable use over extended periods. Adjustable chrome-plated binaurals. Latex-free.	The Uniscope is a single-patient-use stethoscope for situations where cross-contamination is a concern.



Welch Allyn Stethoscopes

Hear the Difference.



4341 State Street Road, PO Box 220, Skaneateles Falls, NY 13153-0220 USA
(p) 1.800.535.6663 (f) 315.685.3361 www.welchallyn.com

WelchAllyn[®]
Advancing Frontline Care™

WelchAllyn[®]
Advancing Frontline Care™

Welch Allyn Stethoscopes

Hear the difference from the leader in diagnostic care

For almost 100 years, we've provided exceptional performance and value from every Welch Allyn device. Our complete lineup of stethoscopes is no exception. From our Harvey DLX with trumpet brass construction for improved sound, to our well-known disposable line, every Welch Allyn product has been designed with you—the clinician—in mind.

Cardiology Grade Stethoscopes

Developed during a 30-year association with Dr. W. Proctor Harvey, MD, a recognized "Master" by the American College of Cardiology.

Welch Allyn Harvey DLX™ Stethoscope

With trumpet brass construction, the Harvey DLX stands out among all other stethoscopes. Trumpet brass construction lets you hear more of the important cardiac, pulmonary, and vascular sounds without relying on difficult "variable pressure" techniques.

- Full-range cardiology/pulmonology stethoscope
- Available in double- or triple-head designs
- Flat or corrugated diaphragms
- Pediatric chestpiece available
- Chrome-plated brass chestpiece and binaurals
- Rotatable binaurals and interchangeable comfort sealing eartips
- Dual-lumen tubing creates individual, uninterrupted sound channels to each ear
- Latex-free
- 10-year warranty

Welch Allyn Harvey™ Elite® Stethoscope

The Harvey Elite offers updated styling with the superior comfort and auscultation capabilities cardiologists have come to rely on. Its stainless steel head transmits sounds crisply and clearly.

- Full-range cardiology/pulmonology stethoscope
- Double head (bell and diaphragm)
- Pediatric chestpiece available
- Stainless-steel chestpiece with non-chill rim
- Rotatable binaurals and interchangeable comfort sealing eartips
- Dual-lumen tubing creates individual, uninterrupted sound channels to each ear
- Flat or corrugated diaphragms
- Latex-free
- 10-year warranty



HARVEY DLX DOUBLE-HEAD STETHOSCOPE 5079-325



HARVEY ELITE DOUBLE-HEAD STETHOSCOPE 5079-125

Welch Allyn Professional Grade Stethoscopes

Excellent acoustic clarity and durability

Providing excellent sound quality, durability and comfort for the practicing professional—they're lightweight and can be custom fitted via adjustable binaurals and eartips—perfect for listening to cardiac, pulmonary, and vascular sounds.

Professional Adult


- Excellent for general auscultation and blood pressure
- Double-head design
- Flat diaphragm
- Adult and Pediatric versions available
- Stainless steel chestpiece and binaurals
- Rotatable binaurals and interchangeable comfort sealing eartips
- Latex-free
- 5-year warranty

Lightweight

- Excellent for general auscultation and blood pressure
- Chrome-plated chestpiece for durability
- Rotatable binaurals
- Double-head design
- Flat diaphragm
- Adult version only
- Latex-free
- 2-year warranty



PROFESSIONAL ADULT STETHOSCOPE 5079-135

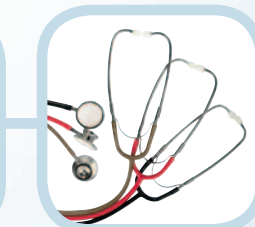

LIGHTWEIGHT STETHOSCOPE 5079-73

Welch Allyn Disposable Stethoscope (Single Patient Use)

Helps protect against cross-contamination

A single-patient-use stethoscope for use in situations where cross-contamination is a concern. Lightweight with excellent acoustics, the Uniscope allows for accurate auscultation without compromising primary stethoscopes or patient safety.

Uniscope

- Single-patient-use (disposable) stethoscope
- Excellent for general auscultation
- Adult (yellow) and Pediatric (blue) versions available
- Plastic chestpiece and binaurals
- Single-head design
- Flat diaphragm
- Latex-free



UNISCOPE 17461